

Rotadorm Care

The classic stand-up bed

Max.
patient
weight
135 kg

Safe
working load
175 kg

In 3 lengths:
200, 210 and
220 cm

Direction
of rotation
possible on
both sides

Height
adjustment
range
40 - 80 cm



◀ With the Rotadorm bed series, the positions of the individual elements of the lying surface can be individually adjusted and programmed to the patient's needs


A stable patient lifting pole and the adjustable ergonomic handle are included in the standard version of our stand-up beds.

The **Rotadorm Care** is the proven care bed with integrated stand-up function from ISKMed. With its high level of functionality, this registered aid provides support and relief in everyday life for both patients and carers.

Product details

- | 85 cm wide single bed, available in 3 lengths
- | Height, head section, neck section, thigh support, lower leg support and turning mechanism individually electronically adjustable
- | 2-section special mattress with incontinence cover
- | Patient lifting pole with handle
- | Removable patient handles
- | Side rails with function lock
- | Programming and patient hand control
- | Stable scissor lift system for height adjustment

Product images

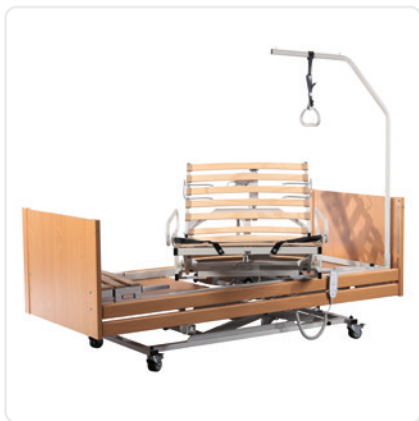

The product images on this page show the standard version of the stand-up bed. The final product changes visually due to optional equipment ("Extra equipment/optional" on next page), such as an additional table. If you have any questions regarding an order, please do not hesitate to contact us .

Technical data

- | Mattress width: 85 cm
- | Mattress length: 200 / 210 / 220 cm
- | Bed height: 40 – 80 cm
- | Max. patient weight: 135 kg
- | Safe working load: 175 kg

Extra equipment / optional

- | Lying surface dimension: 85 x 200 cm (SI-011-1)
- | Lying surface dimension: 85 x 210 cm (SI-011-21)
- | Lying surface dimension: 85 x 220 cm (SI-011-22)
- | Turningside: left
- | Turningside: right
- | **Pull handle (SI-011.85.307)**
serves as a positioning aid in the bed, is mounted on the headboard
- | **Cover for handle receptacle (SI-011.80.725)**
enables getting out of bed to the side
- | **Table (SI-132-0)**
serves as storage option during sitting posture, mounting on grab rail over wire frame, dimension 365 x 265 x 10 mm
- | **Transport rack (TRANSP-08)**
to assist with transport, deposit € 300
- | **Battery pack with charging circuit (SI-011.85.050)**
safety equipment for power failure – all electrical functions can thus be returned to the initial position even without mains voltage