

Rotadorm Low

The stand-up bed for smaller patients



Max.
patient
weight
135 kg

Safe
working load
175 kg

3 lengths
adjustable:
200, 210 and
220 cm

Height
adjustment
range
34 - 78 cm

Direction
of rotation
possible on
both sides



◀ The mattress base can be turned in and out manually using the emergency release on the rotary motor.

▶ The special low-floor design enables particularly low entry and exit and is therefore particularly suitable for smaller patients under 165 cm tall.





The **Rotadorm Low** was developed as a stand-up bed for people under 165 cm tall. It offers all the functionalities of the Rotadorm series, such as individually adjustable frame elements and the option of saving all individual settings for daily use.


Product details

- | 85 cm wide single bed, adjustable to three lengths
- | For individual stand-up procedure with low seat height
- | Height, head section, neck section, thigh support, lower leg support and turning mechanism individually electronically adjustable
- | 2-section special mattress with incontinence cover
- | Patient lifting pole with handle
- | Removable patient handles in frame colour
- | Programming and patient hand control
- | Optionally with battery system for self-sufficient use of all functions
- | Stable scissor lift system for height adjustment
- | Individually braked double swivel castors



Product images

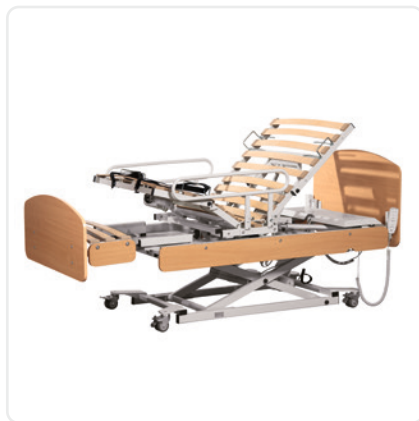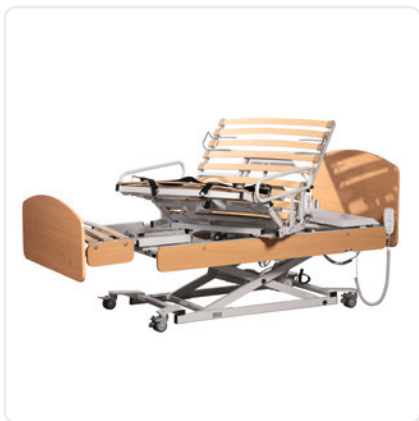

The product images on this page show the standard version of the stand-up bed. The final product changes visually due to optional equipment ("Extra equipment/optional" on next page), such as a pull handle. If you have any questions regarding an order, please do not hesitate to contact us .

Technical data

- | Mattress width: 85 cm
- | Mattress length: 200 / 210 / 220 cm
- | Bed height: 34 – 78 cm
- | Max. patient weight: 135 kg
- | Safe working load: 175 kg

Extra equipment / optional

- | Lying surface dimension: 85 x 200 cm (SI-012-1)
- | Lying surface dimension: 85 x 210 cm (SI-012-21)
- | Lying surface dimension: 85 x 220 cm (SI-012-22)
- | Turningside: left
- | Turningside: right
- | **Pull handle (SI-011.85.307)**
serves as a positioning aid in the bed, is mounted on the headboard
- | **Cover for handle receptacle (SI-011.80.725)**
enables getting out of bed to the side
- | **Table (SI-132-0)**
serves as storage option during sitting posture, mounting on grab rail over wire frame, dimension 365 x 265 x 10 mm