

Exaneuro


The highly stable treatment couch

Large
one-piece
lying surface

Optionally
with **hydraulic**
height
adjustment

Lifting
columns with
electric height
adjustment

Individually
configurable

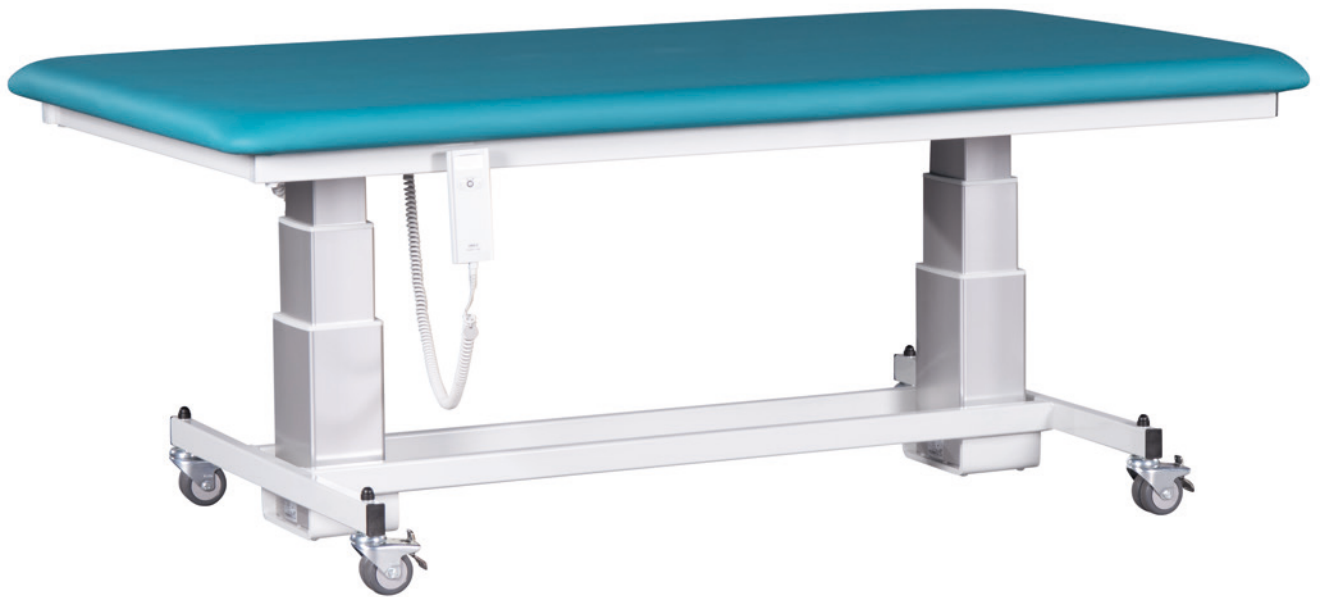

◀ *The column lifting system is driven by powerful, high-quality motors.*

▶ *The columns have a height adjustment range of 40 cm and therefore cover a wide range for physiotherapy treatment with high stability.*




The **Exaneuro** is primarily designed for use with forms of therapy that require a large contact surface and increased stability, such as Vojta or Bobath. The couch can be customised in terms of lying surface dimensions, lifting system and design. All variants are equipped with an integrated safety system.

Product details

- | Lying surface dimensions freely selectable
- | Lifting column
 - Width: 60 – 140 cm
 - Length: 120 – 220 cm
 - Height adjustment range: 48 – 88 cm
 - Under-couch clearance: 10 cm
- | Lifting lever electric
 - Width: 100 – 140 cm
 - Length: 200 – 220 cm
 - Height adjustment range: 44 – 82 cm
 - Under-couch clearance: 13 cm
- | Scissors hydraulic or electric
 - Width: 100 – 120 cm
 - Length: 160 – 220 cm
 - Height adjustment range: 46 – 84 cm
 - Under-couch clearance: 13 cm
- | For all models: height increase possible, please enquire
- | Safe working load up to 300 kg
- | Lifting and propulsion system based on the principle of integrated safety
- | Single braked swivel castors

Product images



The first product image series on this page shows the treatment couch design in the lifting column version. After the color variants, you will also see the couch in the “scissors electric” version and subsequently “lifting lever electric”. The final product consequently changes optically due to optional equipment (“Extra equipment/optional” on next page), such as a roller lifting system. If you have any questions regarding an order, please do not hesitate to contact us [contact us](#) .

Note on colour representation: The product images shown on this page are for illustrative purposes only. Colour variations between the screen display and the actual product are possible due to different display technologies and individual device settings. This applies in particular to the metal colour RAL 7035 shown. Only the original RAL colour tone according to the manufacturer’s specifications is authoritative.

Extra equipment / optional

| Lifting column



| Lifting lever electric



| Scissors hydraulic



| Scissors electric



| Additional foot switch for adjusting the working height (BL-101-2)

| Nose slit in the lying surface or head section (SB-187-0)

| Cover for nose slit (SB-187-1)

| Head section electrically adjustable; +/- 30° (BL-312-0)
(surcharge for ratchet adjustment)

| Side rails for treatment couch, chromed (1 pair) (BL-201-0)

| Side rails for treatment couch, chromed (1 piece) (BL-202-0)

| Head section adjustable with ratchet assembly; +/- 30° (BL-310-0)

| Head section adjustable with gas spring (BL-311-0)
(surcharge for ratchet adjustment)

| Leg section with ratchet adjustment (BL-999-3)
Foot section mechanically adjustable up to +30° with ratchet adjuster

| Paper roll holder (BL-151-0)

| Roller lifting system (BL-501-0)

| Central brake system

| Padded/foldable surround - for up to 4 sides

| Pacificblue (No. 432x1385)
Colour imitation leather cover



| Mediterranean (No. 208x1752)
Colour imitation leather cover



| Saffron (No. 208x2388)
Colour imitation leather cover



| Racingred (No. 207x1173)
Colour imitation leather cover



| Mercedes Date (No. 207x3115)
Colour imitation leather cover



| Special colour
Example: Skai Tundra
Wide range of additional colours

