

Exatrain Flex

The flexible stand-up board

Individually
configurable

Electric
height
adjustment:
44 - 84 cm

Swivelling
range
up to 85°

Belt
support
system
included



◀ In the event of a power failure, the couch can be returned to the horizontal position manually at any time by activating the emergency release.

▶ The belt support system holds the patient securely in the standing position.



The **Exatrain Flex** is a special couch for standing training, which is used in physiotherapy, neurology, surgery and orthopaedics, among others. The highly stable substructure and the belt support system provide the patient with absolute safety during training in any position. In addition to the standard version, the Exatrain Flex can be configured with a wide range of accessories according to the needs and requirements of the practitioner.

Product details




- | Lying surface dimensions: 73 x 185 cm
- | Under-couch clearance: 13 cm (height increase possible, please inquire)
- | Maximal patient weight: 135 kg
- | Safe working load: 170 kg
- | Footboard fully leathered
- | Optionally with adjustable electric head section
- | Various fixing and table systems available as accessories
- | Individually braked double castors
- | Integrated safety through automatic locking of electrical functions



Product enquiry
Exatrain Flex

**| Aid number for German
health insurance companies:
28.29.02.1046**

Product images

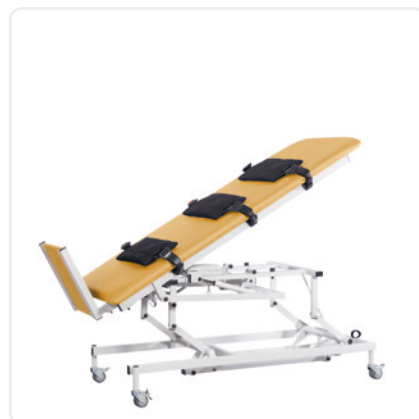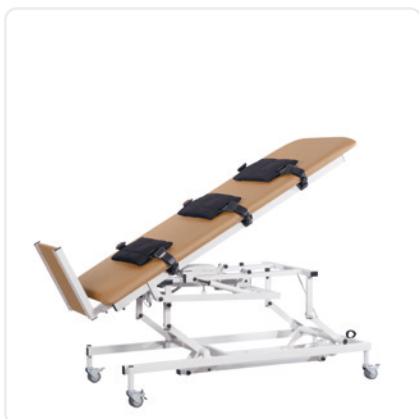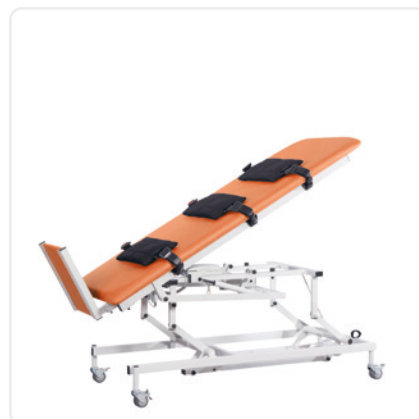

The product images on this page show the standard version of the standing board. The final product changes visually due to optional equipment ("Extra equipment/optional" on next page), such as metal side rails. If you have any questions regarding an order, please do not hesitate to **contact us** .

Note on colour representation: The product images shown on this page are for illustrative purposes only. Colour variations between the screen display and the actual product are possible due to different display technologies and individual device settings. This applies in particular to the metal colour RAL 7035 shown. Only the original RAL colour tone according to the manufacturer's specifications is authoritative.

Technical data

- | Width: 73 – 120 cm
- | Length: 185 – 205 cm
- | Height: 44 – 84 cm
- | Safe working load: 170 kg

Extra equipment / optional

- | Footplates individually angle and length adjustable (SB-121-0)
- | Metal side rails
- | Armrests (SB-131-0)
- | Working and therapy table with body cut-out (wood), 90 cm (SB-141-0)
- | Working and therapy table with body cut-out (plexiglass), 90 cm (SB-141-1)
- | Electrical head section adjustment (SB-170-0)
- | Slit in the lying surface for single fixation of the legs (SB-171-0)
- | Head fixation with mounting for lying surfaces up to 100 cm width (SB-181-0)
- | Protractor (angle gauge) for mounting on the lying surface (SB-185-0)
- | Nose slit for stand-up board (8 cm x 18 cm) (SB-187-0)
- | Cover for nose slit (8 cm x 18 cm) (SB-187-1)
- | Nose slit for stand-up board (8 cm x 30 cm) (SB-187-2)
- | Cover for lying surface length slit (SB-187-4)
- | Chassis elevation for underpassing with lifter (SB-188-0)
- | Imitation leather belt system (SB-151-5)
- | Battery system



Product enquiry
Exatrain Flex

ISKO Koch GmbH · Egerländer Str. 28 · 95448 Bayreuth
Tel +49 921 150845-0 · info@iskomed.de · www.iskomed.com